



Test Report

Date : 2020-07-23
No. : HC20070472

Page 1 of 6

Applicant(Code:01325270) : CareHK Limited
Unit 2A Precision Manufacturing Centre
3 Dai Hei Street
Tai Po Industrial Estate NT HK


Description of Sample(s) : One submitted sample said to be CareHK Face Mask.
Item No.: CHK001-AH / CHK001-KH
Batch / Lot No.: A191011
Country of Origin & Destination : Hong Kong

Sample(s) Received Condition(s): In intact original package under ambient temperature

Date Sample(s) Received : 2020-07-13

Date Tested : 2020-07-16 to 2020-07-23

Investigation Requested : Performance Test as per ASTM F2100
1. Bacterial Filtration Efficiency 細菌過濾效率 (BFE) %
– *Staphylococcus aureus* 金黃色葡萄球菌 (ATCC 6538)
2. Particulate Filtration Efficiency 微粒過濾效率 (PFE) %
3. Differential Pressure 空氣交換壓力差
4. Synthetic Blood Penetration 合成血液穿透性
5. Flammability to Class 1 可燃性


CHUAH Ah Poh, Anne
Authorized Signatory
Chemical and Food Department
For and on behalf of
The Hong Kong Standards and Testing Centre Ltd.



In the event of any discrepancy between the English and Chinese versions of the text mentioned before, the English version shall prevail to the extent of such discrepancy only.

The Hong Kong Standards and Testing Centre Limited 香港標準及檢定中心有限公司
10 Dai Wang Street, Tai Po Industrial Estate, Tai Po, N.T., Hong Kong 香港新界大埔工業邨大宏街10號
Tel: +852 2666 1888 Fax: +852 2664 4353 Email: hkstc@stc.group Website: www.stc.group

This report shall not be reproduced unless with prior written approval from The Hong Kong Standards and Testing Centre Limited.

For Conditions of Issuance of this test report, please refer to the overleaf and Website.

除非事先得到香港標準及檢定中心有限公司的書面批准，否則不得自行複印此報告。測試報告相關的簽發條款，請參閱背頁和本公司網頁。



Test Report

Date : 2020-07-23
No. : HC20070472

Page 2 of 6

Requirement:

Performance Test as per ASTM F2100-19	Level 1	Level 2	Level 3
Bacterial Filtration Efficiency(BFE) % 細菌過濾效率 – <i>Staphylococcus aureus</i> (ATCC 6538) 金黃色葡萄球菌	≥95%	≥98%	
Particulate Filtration Efficiency (PFE) % 微粒過濾效率	≥95%	≥98%	
Differential Pressure (ΔP) 空氣交換壓力差	<5.0 mmH ₂ O/cm ²	<6.0 mmH ₂ O/cm ²	
Resistance to Penetration by Synthetic Blood Penetration 合成血液穿透性	80 mmHg	120 mmHg	160 mmHg
Flame Speed (Flammability to Class 1) 可燃性	Class 1 (The time of flame spread is 3.5 seconds or more)		

Summary:

Performance Test as per ASTM F2100-19	CareHK Face Mask Batch / Lot No.: A191011
	Level 3
Bacterial Filtration Efficiency(BFE) % 細菌過濾效率 – <i>Staphylococcus aureus</i> (ATCC 6538) 金黃色葡萄球菌	Pass
Particulate Filtration Efficiency (PFE) % 微粒過濾效率	Pass
Differential Pressure (ΔP) 空氣交換壓力差	Pass
Resistance to Penetration by Synthetic Blood Penetration 合成血液穿透性	Pass
Flame Speed (Flammability to Class 1) 可燃性	Pass

The Hong Kong Standards and Testing Centre Limited 香港標準及檢定中心有限公司
10 Dai Wang Street, Tai Po Industrial Estate, Tai Po, N.T., Hong Kong 香港新界大埔工業邨大宏街10號
Tel: +852 2666 1888 Fax: +852 2664 4353 Email: hkstc@stc.group Website: www.stc.group

This report shall not be reproduced unless with prior written approval from The Hong Kong Standards and Testing Centre Limited.

For Conditions of Issuance of this test report, please refer to the overleaf and Website.

除非事先得到香港標準及檢定中心有限公司的書面批准，否則不得自行複印此報告。測試報告相關的簽發條款，請參閱背頁和本公司網頁。



Test Report

Date : 2020-07-23
No. : HC20070472

Page 3 of 6

Method(s) Used:

1. ASTM F2100-2019 9.1 & ASTM F2101-2019
2. ASTM F2100-2019 9.3 & ASTM F2299-2017
3. ASTM F2100-2019 9.2 & EN 14683:2019 Annex C
4. ASTM F2100-2019 9.4 & ASTM F1862-2017
5. 16 CFR 1610

Test Result(s):

1. Bacterial Filtration Efficiency 細菌過濾效率 (BFE) %
– *Staphylococcus aureus* 金黃色葡萄球菌 (ATCC 6538)

Specimen(s)	CareHK Face Mask Batch / Lot No.: A191011
1	99.3%
2	99.3%
3	99.1%
4	99.2%
5	99.0%

Notes : - Control average 對照組的平均菌落數 : 1842CFU
- Mean particle size 平均粒徑 : 2.8µm
- Testing side 測試面 : Outside of specimen 外側
- Testing area 測試面積 : 47 cm²

The Hong Kong Standards and Testing Centre Limited 香港標準及檢定中心有限公司
10 Dai Wang Street, Tai Po Industrial Estate, Tai Po, N.T., Hong Kong 香港新界大埔工業邨大宏街10號
Tel: +852 2666 1888 Fax: +852 2664 4353 Email: hkstc@stc.group Website: www.stc.group

This report shall not be reproduced unless with prior written approval from The Hong Kong Standards and Testing Centre Limited.

For Conditions of Issuance of this test report, please refer to the overleaf and Website.

除非事先得到香港標準及檢定中心有限公司的書面批准，否則不得自行複印此報告。測試報告相關的簽發條款，請參閱背頁和本公司網頁。



Test Report

Date : 2020-07-23
No. : HC20070472

Page 4 of 6

2. Particulate Filtration Efficiency 微粒過濾效率 (PFE) %
(0.1µm PSL ; Flow rate : 28 Litre/min)

Specimen(s)	CareHK Face Mask Batch / Lot No.: A191011
1	99.6%
2	99.8%
3	99.6%
4	99.4%
5	99.7%
Average	99.6%

3. Differential Pressure 空氣交換壓力差

Specimen(s)	CareHK Face Mask Batch / Lot No.: A191011
1	3.9 mmH ₂ O/cm ²
2	4.3 mmH ₂ O/cm ²
3	4.5 mmH ₂ O/cm ²
4	4.6 mmH ₂ O/cm ²
5	4.2 mmH ₂ O/cm ²

Notes : - 1mmH₂O = 9.8Pa

The Hong Kong Standards and Testing Centre Limited 香港標準及檢定中心有限公司
10 Dai Wang Street, Tai Po Industrial Estate, Tai Po, N.T., Hong Kong 香港新界大埔工業邨大宏街10號
Tel: +852 2666 1888 Fax: +852 2664 4353 Email: hkstc@stc.group Website: www.stc.group

This report shall not be reproduced unless with prior written approval from The Hong Kong Standards and Testing Centre Limited.

For Conditions of Issuance of this test report, please refer to the overleaf and Website.

除非事先得到香港標準及檢定中心有限公司的書面批准，否則不得自行複印此報告。測試報告相關的簽發條款，請參閱背頁和本公司網頁。



Test Report

Date : 2020-07-23
No. : HC20070472

Page 5 of 6

4. Synthetic Blood Penetration 合成血液穿透性 (Pressure : 160mmHg)

Specimen(s)	CareHK Face Mask Batch / Lot No.: A191011	Specimen(s)	CareHK Face Mask Batch / Lot No.: A191011
1	None Seen 無穿透	17	None Seen 無穿透
2	None Seen 無穿透	18	None Seen 無穿透
3	None Seen 無穿透	19	None Seen 無穿透
4	None Seen 無穿透	20	None Seen 無穿透
5	None Seen 無穿透	21	None Seen 無穿透
6	None Seen 無穿透	22	Penetration 有穿透
7	None Seen 無穿透	23	None Seen 無穿透
8	None Seen 無穿透	24	None Seen 無穿透
9	None Seen 無穿透	25	None Seen 無穿透
10	Penetration 有穿透	26	None Seen 無穿透
11	None Seen 無穿透	27	None Seen 無穿透
12	None Seen 無穿透	28	None Seen 無穿透
13	None Seen 無穿透	29	None Seen 無穿透
14	None Seen 無穿透	30	None Seen 無穿透
15	None Seen 無穿透	31	None Seen 無穿透
16	None Seen 無穿透	32	None Seen 無穿透
Requirement: An acceptable quality limit of 4.0% is met for a normal single sampling plan when ≥ 29 of 32 test specimens show passing result (none seen)			

The Hong Kong Standards and Testing Centre Limited 香港標準及檢定中心有限公司
10 Dai Wang Street, Tai Po Industrial Estate, Tai Po, N.T., Hong Kong 香港新界大埔工業邨大宏街10號
Tel: +852 2666 1888 Fax: +852 2664 4353 Email: hkstc@stc.group Website: www.stc.group

This report shall not be reproduced unless with prior written approval from The Hong Kong Standards and Testing Centre Limited.

For Conditions of Issuance of this test report, please refer to the overleaf and Website.

除非事先得到香港標準及檢定中心有限公司的書面批准，否則不得自行複印此報告。測試報告相關的簽發條款，請參閱背頁和本公司網頁。



Test Report


Date : 2020-07-23
No. : HC20070472

Page 6 of 6

5. Flammability

Specimen(s)	CareHK Face Mask Batch / Lot No.: A191011	Class
	Time of spread of flame (Original state)	
1	Did not ignite 沒有燃點	1
2	Did not ignite 沒有燃點	
3	Did not ignite 沒有燃點	
4	Did not ignite 沒有燃點	
5	Did not ignite 沒有燃點	

Photo(s):

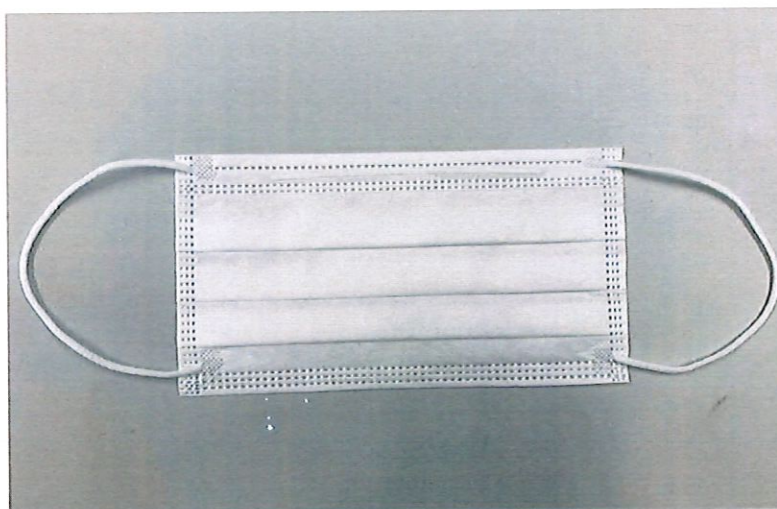

***** End of Test Report *****

The Hong Kong Standards and Testing Centre Limited 香港標準及檢定中心有限公司
10 Dai Wang Street, Tai Po Industrial Estate, Tai Po, N.T., Hong Kong 香港新界大埔工業邨大宏街10號
Tel: +852 2666 1888 Fax: +852 2664 4353 Email: hkstc@stc.group Website: www.stc.group

This report shall not be reproduced unless with prior written approval from The Hong Kong Standards and Testing Centre Limited.

For Conditions of Issuance of this test report, please refer to the overleaf and Website.

除非事先得到香港標準及檢定中心有限公司的書面批准，否則不得自行複印此報告。測試報告相關的簽發條款，請參閱背頁和本公司網頁。